


DiaSpect Tm

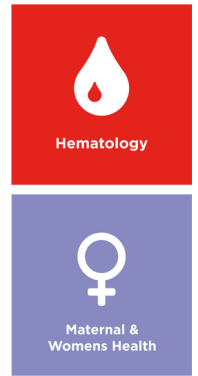

DiaSpect Tm

Immediate hemoglobin results
from a hand-held analyzer

- ◆ Results in just about 1 second
- ◆ Reagent-free microcuvettes
- ◆ 40 days or 10,000 tests continuous use
- ◆ 2½ years cuvette shelf life

Fast, easy to use and reliable in all climate conditions

<p>Fast and accurate</p>	<ul style="list-style-type: none"> • Results in just about 1 second • Measuring range: 0 - 25.5 g/dL • Imprecision (within run): CV < 1% • Factory calibrated against HiCN reference method according to ICSH • Photometric method with broad spectrum and compensation for turbidity and scattering • No additional calibration required
<p>Easy to use</p>	<ul style="list-style-type: none"> • 'Always on' technology. The analyzer is always ready for use • Automatic self-check between every measurement • Long microcuvette shelf life: 2.5 years from date of production • Microcuvettes unaffected by humidity or temperature and can be used until expiration date after opening the bag
<p>Practical and portable</p>	<ul style="list-style-type: none"> • Sample volume: <10 μL • Capillary, venous or arterial blood (EDTA and Heparin anticoagulant can be used) • Maintenance free • Operating ambient temperature: 10 °C to 42 °C • Storage temperature: 0 °C to 50 °C (- 30 °C to 70 °C up to 24 hours maximum.) • Integrated rechargeable lithium-ion battery (40 days / 10,000 tests of continuous use) • Compact: 150 (d) x 90 (w) x 40mm (h) • Lightweight: 180 g • Connect to PC via USB 2.0 or Bluetooth (optional)



1

Collect blood sample



2

Insert the filled cuvette and press down gently



3

Results appear in just about 1 second

Distributed by

Manufacturer
 DiaSpect Medical GmbH
 Ebendorfer Chaussee 3
 39179 Barleben
 Germany
 +49 (0) 39203 511 0
 info@ekf-diagnostic.de
 www.ekfdiagnostics.com

EKF | Diagnostics
 for life